

## The Sonesta line of tables and chairs



Product catalog for models

6202 | 6303 | 6302 | 6302U | 6210

# Gynecology/Urology Procedure Chair

## SONESTA 6202



### SPECIFICATIONS:

LENGTH WITH EXTENSION CUSHION	190cm	75"
LENGTH WITHOUT	143cm	56"
WIDTH	62cm	25"
HEIGHT WITH SEAT CUSHION (MIN)	40cm	16"
HEIGHT WITH SEAT CUSHION (MAX)	95cm	37"
LIFT CAPACITY	170kg	375lbs
WHEEL DIAMETER (LOCKING)	10cm	4"
NET WEIGHT	<125kg	<276lbs

### ELECTRICAL:

120/230v 24v DC Output  
Meets EMC (Electromagnetic Compatibility) for IEC 60601  
All Products carry the CE mark.

### The ultimate in positioning for your practice.


The Sonesta 6202 is the ultimate chair for Urodynamics and Any Urology examination. Popular accessories, such as fold up foot rests and a seat cushion without cut out are included as standard. Additionally the locking mechanism for one touch movement, offers best-in-class comfort for the patient and the best available ergonomics and ease-of-use for the physician.

### Engineered to be the best, most flexible solution in your examination room:

- Locking mechanism for tensionless and precise positioning
- Three motors for height, seat and back adjustments
- Easy positioning for elderly or physically challenged patients
- Wide range of adjustment, from a low of 16" (40cm) to 39.4" (100cm) high
- Engineered ergonomics. Infinitely adjustable (see specifications) for physician and patient comfort.
- Lightweight and mobile easy to move to other rooms; unmatched access for floor cleaning and disinfecting.
- Tough and reliable, able to support up to 375 lbs-170kg for years of service in the most demanding practices.
- Choose between several color schemes for a perfect fit in your office.

### Standard Accessories:

3 motor control/foot plate	524-6358
1 cushion with 170mm cutout	524-6366-3
1 cushion with without cutout	524-6365
1 Neck pillow	524-6346
2 fold up footrests	524-6331
2 leg rest holders	524-6338
2 heavy duty Rissler leg rests	524-6336
2 hand grips/arm rests	524-6339-2
Drainage receptacle	524-6325-1
Paper roll holder	524-6326

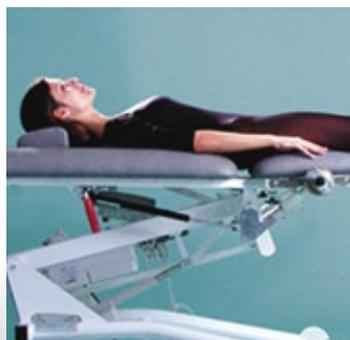

# Fluoroscopy Procedure Table

## SONESTA 6210



### SPECIFICATIONS:

LENGTH WITHOUT FOOT REST/LEG EXTENSION	156cm	63"
LENGTH WITH FOOT REST/LEG EXTENSION	196cm	78"
WIDTH (TOTAL):	99cm	39"
HEIGHT WITH SEAT CUSHION (MIN):	50cm	20"
HEIGHT WITH SEAT CUSHION (MAX):	130cm	51"
IMAGING AREA WITHOUT FOOT REST/LEG EXTENSION (LENGTH)	156CM	63"
IMAGING AREA WITH FOOT REST/LEG EXTENSION (LENGTH)	196CM	78"
WEIGHT CAPACITY	250KG	551LBS

### ELECTRICAL:

120/230v 24v DC Output  
Meets EMC (Electromagnetic Compatibility) for IEC 60601  
All Products carry the CE mark.

### A versatile examination table with C-arm capabilities.

The 6210 fluoroscopy procedure table offers the physician standing, seated, and supine positioning capabilities. This versatile table easily adjusts from/to standing, seated, and supine positions via a hand held control, thus creating a table that is easily accessible to patients of any age.

The side-mounted cantilevered design offers the physician anterior/posterior, lateral and oblique views without obstruction. For added convenience the table is equipped with a battery backup for short-term use in case of power failures or in emergency situations. The superb c-arm access and a wide range of positioning possibilities, make the 6210 the ideal table for any fluoroscopy procedure.



### Positioned to be the best:

- Three motors offer unparalleled positioning/ access options
- Radiolucent carbon fiber seat and back for head-to-toe imaging
- 4 Double wheel rolling casters offer unencumbered maneuverability

### Standard Accessories:

Hand control	524-6356
1 plastic funnel w/drainage	6325-9
1 Fixed arm rest	524-6021-6

# Gynecology/Urology Procedure Chair/Table

SONESTA 6302

SONESTA 6302U

SONESTA 6303



## SPECIFICATIONS:

LENGTH WITH EXTENSION CUSHION	175cm	69"
LENGTH WITHOUT	125cm	49.2"
WIDTH	65cm	25.5"
HEIGHT WITH SEAT CUSHION (MIN)	40cm	16"
HEIGHT WITH SEAT CUSHION (MAX)	100cm	39"
LIFT CAPACITY	170kg	375lbs
WHEEL DIAMETER (LOCKING)	10cm	4"
NET WEIGHT	<100kg	<220lbs

## ELECTRICAL:

120/230v 24v DC Output


Meets EMC (Electromagnetic Compatibility) for IEC 60601

All Products carry the CE mark.

## A flexible solution in your examination room.

The Sonesta 6300 series is a flexible solution to your examination room. With three motor functionality the model 6303 will easily adjust the position of your patient with lift, reclining back and a pelvic tilt, all at the touch of a button. Two motor maneuvering is available with an adjustable back and lift motor in model 6302U, where the seat portion is fixed to one of three desired positions. Model 6302 incorporates motors into the lift and seat section, and uses hydraulic pistons to assist in manually raising the back portion of the procedure chair. The Sonesta 6300 series is packaged with basic accessories yet is compatible with the full line of Sonesta add-ons making for a flexible, and comfortable procedure environment.



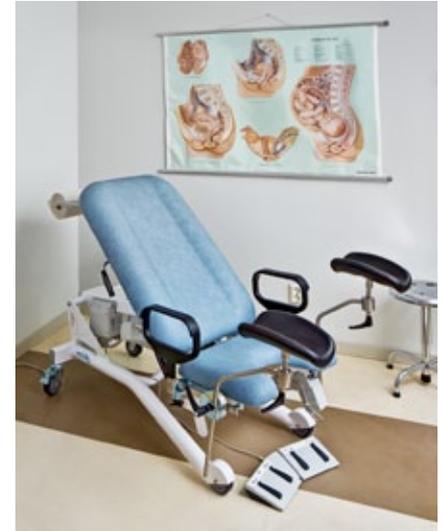

## SONESTA 6303

### Your Ergonomic chair with Procedure versatility:

The 6303 is an affordable path to the unique Sonesta flexibility. Three motors provide full height, seat, and back adjustments. Packaged with basic accessories, the 6303 accepts the full line of Sonesta add-ons, allowing you to create the perfect table for your needs.

The 6303's ergonomic design delivers a higher level of patient and physician comfort and flexibility making it an incredible value in today's market!

- Sonesta flexibility for office, clinic, or hospital
- Three Motors for height, seat and back adjustments
- Easy positioning for patient and physician comfort
- Accepts the full range of Sonesta accessories
- Wide range of adjustment, from a low of 40cm (16") to 100cm (39.4") high
- Lightweight and mobile easy to move to other rooms; unmatched access for floor cleaning and disinfecting.
- Tough and reliable supports up to 375 lbs - 170 kg for years of service in the most demanding practices.

### Standard Accessories:

3 motor foot control	524-6357
Cushion w/170mm cutout	524-6366-3
2 leg rest holders	524-6338
2 heavy duty Rissler leg rests	524-6336
2 hand grips/arm rests	524-6339-2
Drainage receptacle	524-6325-1
Paper roll holder	524-6326

## SONESTA 6302

## SONESTA 6302U

### Flexibility for the office on a budget:

This two-motor table offers complete flexibility for offices on a budget. Motorized height & seat position, with gas piston assisted control of back angle. Or, choose the 6302U: fixed seat with three tilt-position settings and two-motored powered height and back control. Both models include a basic accessory package to cover many procedures, and they accept the full line of Sonesta accessories.

- Sonesta flexibility for office, clinic, or hospital
- Easy positioning for patient and physician comfort
- Accepts the full range of Sonesta accessories
- Wide range of adjustment, from a low of 40cm (16") to 100cm (39.4") high
- Lightweight and mobile easy to move to other rooms; unmatched access for floor cleaning and disinfecting.
- Tough and reliable supports up to 375 lbs - 170 kg for years of service in the most demanding practices.

### Standard Accessories:

2 motor foot control	524- 6353 / 524-P6409
Cushion w/170mm cutout	524-6366-3
2 leg rest holders	524-6338
2 heavy duty Rissler leg rests	524-6336
2 hand grips/arm rests	524-6339-2
Drainage receptacle	524-6325-1
Paper roll holder	524-6326

# Features Comparison Chart



## 6302 your Basic chair for urodynamics

The Sonesta 6302/6302U meets your needs for daily urology and gynecology exams.

## 6303 your Classic chair for urodynamics

With 3 motor controls the Sonesta 6303 will give you all the positions necessary to see a wide-variety of patients.



## 6202 your Ultimate chair for urodynamics

With 3 motor controls along with a locking mechanism, the Sonesta 6202 will give you all the positions necessary to see a wide-variety of patients.




## 6210 Engineered to Change Your Point-of-View

The 6210 video fluoroscopy procedure table is a perfect economic alternative for all your c-arm procedures. With superior diagnostic accuracy, and a variety of patient positioning options, it is a perfect fit for any practice.

	6202	6303	6302	6302U	6210
Feature	Gynecology/Urology Procedure Chair	Gynecology/Urology Procedure Chair	Gynecology/Urology Procedure Chair	Gynecology/Urology Procedure Chair	Fluoroscopy Procedure Chair
Easy Positioning for Elderly or Disables	X	X	X	X	X
Motorized Height Adjustment	X	X	X	X	X
Motorized Seat Adjustment	X	X	X		X
Motorized Reclining Back	X	X		X	X
Locking Device (Seat and Back)	X				
Simultaneous adjustable Seat and Back	X				X
Radiolucent Carbon Fiber Seat and Back					X
Patient Capacity 375 lbs - 170 kg	X	X	X	X	
Patient Capacity 551 lbs - 250 kg					X

# Standard Accessories Comparison Chart

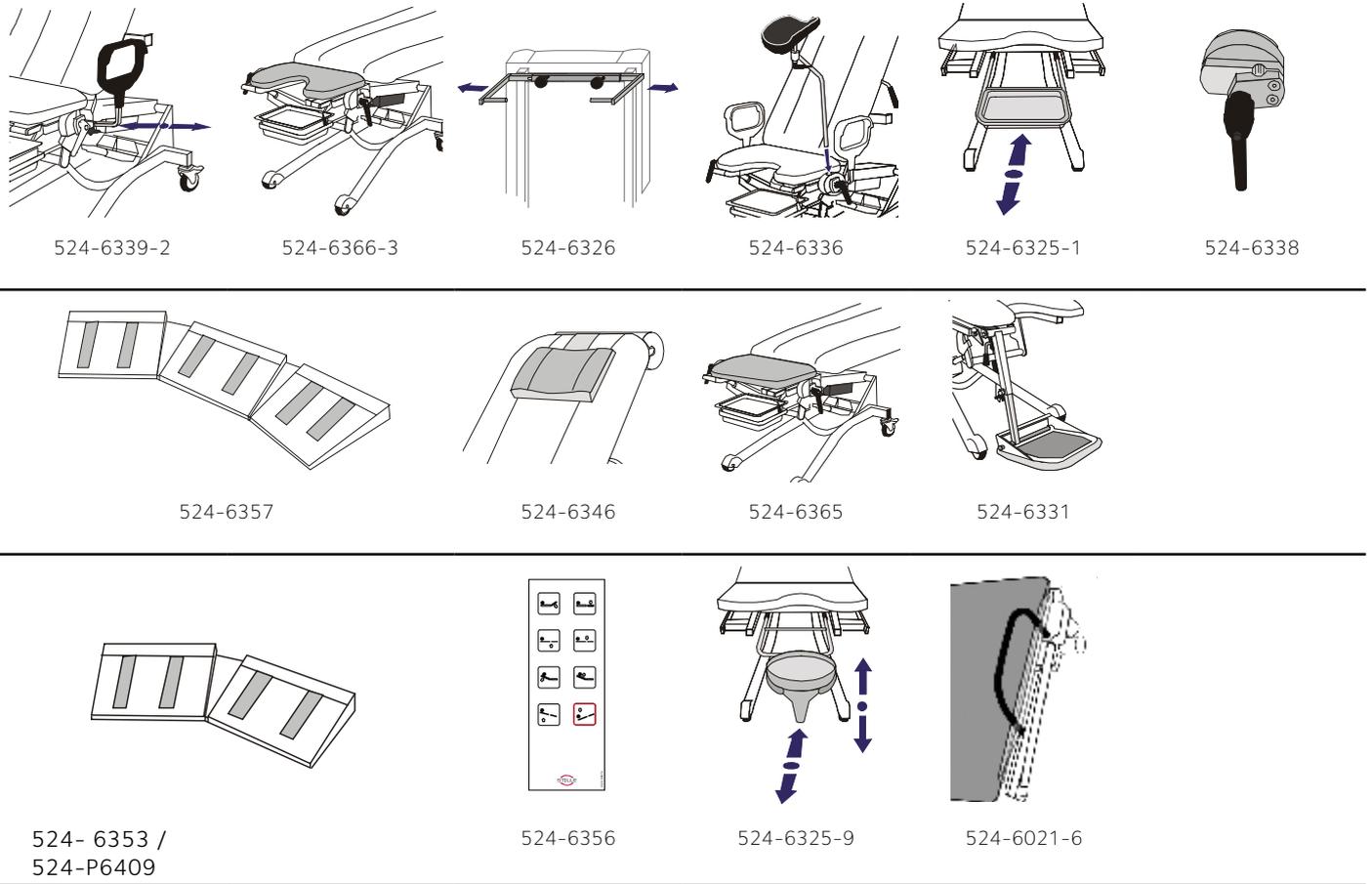

For the complete accessory line, please see the accessory brochure

Part Number	Description	6202	6303	6302	6302U	6210
		Gynecology/Urology Procedure Chair	Gynecology/Urology Procedure Chair	Gynecology/Urology Procedure Chair	Gynecology/Urology Procedure Chair	Fluoroscopy Procedure Chair
524-6339-2	Patient Arm Rests	X	X	X	X	
524-6366-3	Seat Cushion w/ 170 mm Cutout	X	X	X	X	
524-6326	Paper Roll Holder	X	X	X	X	
524-6336	Heavy Duty Rissler Leg Rests	X	X	X	X	
524-6325-1	Drainage Receptacle	X	X	X	X	
524-6338	Leg Rest Holder	X	X	X	X	
524-6357	Three Motor Foot Control	X	X			
524-6346	Neck Pillow	X				
524-6365	Seat Cushion w/o Cutout	X				
524-6331	Fold Up Footrests	X				X
524-6353	Two Motor Foot control (6302)			X		
524-P609	Two Motor Foot control (6302U)				X	
524-6356	Hand Control for 6210					X
524-6325-9	Plastic Funnel /w Drainage for 6210					X
524-6021-6	Arm Rest for 6210 (1 Fixed)					X



Distributed by:

**Sonesta Medical**

PHONE:  
FAX:  
Email: [info@sonestamedical.com](mailto:info@sonestamedical.com)  
[www.stille.se](http://www.stille.se)